



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Product Manual of DAQIDA

Choking Rescue Device set –Manual Plunger

Suction System for Adults and Children

B2B Product Specification & Description Manual | DAQIDA

Product Overview

The Choking Rescue Device – Manual Plunger Suction System (MODEL: PZHXQ-F01) is a pure mechanical emergency rescue tool designed for wholesale and global B2B distribution. Requiring no batteries or electricity, this portable bottle-style anti-choking device provides immediate airway obstruction relief for adults and children, making it ideal for household, vehicle, and public safety preparedness.



Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Key Product Features

Vent Valve: Releases internal pressure control and works alongside the system to facilitate one-way air flow.

Handle: Ergonomically shaped for a secure grip during emergency operation.

Plunger & Bellows: Designed as a flexible PE corrugated pump body that can be repeatedly compressed and pulled without cracking, generating stable negative pressure.

Face Mask: Ultra-soft, skin-friendly medical PVC material designed to comfortably fit different facial contours while maintaining an airtight seal. Includes interchangeable options for adults and children.



How to Use

Step 1: PLACE: Place the mask over the victim's nose and mouth. Hold the mask with your hand so that it forms a good all-around seal on the victim's face. If the victim is unconscious, carefully place them on a flat and sturdy surface, and choose a mask size suitable for the victim.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Step 2: PUSH: Push the plunger down with your strong hand; all air will escape out of the one-way valve system, preventing an obstruction from being pushed further into the victim. **Step 3:** PULL: Keep the mask sealed and pull upward to generate powerful negative pressure that extracts the foreign object.

Step 3 PULL: Pull back on the plunger while holding the mask in place to create suction. If the obstruction does not come out on the first pull, continue to push and pull the plunger without removing the device from the victim's face until the obstruction is removed. cleared.

MODEL:PZHXQ-F01



*If the victim is unconscious, carefully place it on a flat and sturdy surface.
*Choose a mask size that is suitable for the victim; If used on a child, please make sure to use children's sizes.
*Place the mask on the mouth and nose, ensuring a tight seal to generate effective suction.

EASY TO USE IN CHOKING EMERGENCIES

PLACE


Place the mask over the victim's nose and mouth. Hold the mask with your hand so that it forms a good all-around seal on the victim's face.

PUSH

Push the plunger down with your strong hand, all air will escape out of the one-way valve system preventing an obstruction from being pushed further into the victim.

PULL

Pull back on the plunger while holding the mask place. This will create suction. If the obstruction does not come out on the first pull, continue to push, and pull the plunger without removing the device from the victim's face until the obstruction is removed.

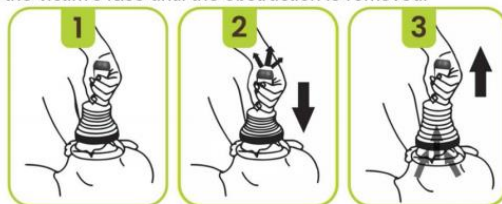

Materials & Technical Selection & Engineering

Face Mask (Soft PVC): Selected for its soft texture to tightly fit various facial shapes, ensuring an effective seal and negative pressure during rescue. The transparent material allows rescuers to monitor the victim's facial condition, while being skin-friendly to reduce irritation.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Bellows, Handle & Base (PE Polyethylene): Features high toughness to withstand repeated push-pull telescoping without cracking, ensuring stable negative pressure generation. Lightweight and stable, it resists aging and embrittlement during long-term emergency storage.

Independent Packaging (PE Film): Dustproof and moisture-proof film that individually protects the mask to maintain cleanliness and ease of carrying and storage.

Design Advantage: Combines PVC and PE to leverage their respective strengths—PVC provides flexible sealing, while PE handles repeated deformation to generate suction.



Application Scenarios & Use Cases

Household Emergency Readiness: An essential safety device for family rooms and kitchens to protect household members from sudden choking incidents.

Automotive & Travel Safety: Compact, portable, and durable design suitable for storage in vehicle glove boxes or travel emergency kits.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Institutional & Public Safety: Reliable, non-invasive safety equipment for schools, daycares, restaurants, and hospitality venues.ar.

SUITABLE FOR A WIDE RANGE OF SCENARIOS

Hands-on training for real choking emergencies.



Easy to Learn



Safe to Use



Quick Response



Life Saving



Technical Specifications

Specification Parameter	Details
Product Name	Choking Rescue Device
Set Specification	Single Pack
Product Color	Transparent green
Product Size (cm)	Adult Mask: 12.5×10.0×7.5CM, Kids Mask: 10×7.5×6CM, Main Body: 16×9.5×9.5CM
Product Package Size (cm)	25×25×10CM

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Product Material	PE+PVC
Product Net Weight (g)	167 g
Product Gross Weight (with Package, g)	179g
Packaging Material	Pc packing
Outer Carton Size (cm)	54×47×47CM
Net Weight per Carton (kg)	12.5 KG
Gross Weight per Carton (kg)	13.55kg
Quantity per Carton (Sets)	70 PCS
Custom Color MOQ (Sets)	1000
Custom Logo MOQ (Sets)	1000
Invoicing	(Please inform in advance) 1% Regular Invoice / 13% VAT Invoice

Patented Choking Rescue Device Set -


United States Design Patent: Protected under US Patent No.: US D1071,383 S (Date of Patent: April 15, 2025), securing the ornamental bottle-style configuration, ergonomic plunger handle, and integrated corrugated bellows geometry.

European Union Intellectual Property Office (EUIPO): Officially registered under Community Design No. 01504989-0001, safeguarding the unique aesthetic and structural form across the European market.

Global Compliance & Safety Standards: Backed by FDA establishment registration, SGS materials testing reports, and rigorous third-party quality certifications ensuring full export readiness and legal protection for brand owners and distributors.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province

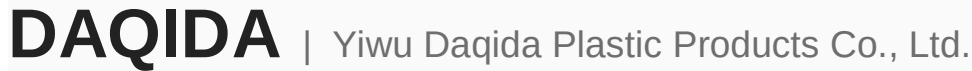

Safety • Quality • Compliance

[illegible]