



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Product Manual of DAQIDA

Choking Rescue Device – Manual Airbag Suction System for Adults and Children

B2B Product Specification & Description Manual | DAQIDA

Product Overview

The Choking Rescue Device – Manual Airbag Suction System for Adults and Children is an OEM-ready emergency clearance kit designed for wholesale and global B2B distribution. As a vital airway clearance device, this portable anti-choking suction pump requires no electricity or batteries, making it an essential wholesale safety product for medical distributors, school safety programs, and household emergency kits.



Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Key Product Features

Structural Stability: Built with a resilient mechanical design for smooth, bind-free linear motion during bulk emergency deployment.

Integrated One-Way Valve: Protects patients by expelling air outward without pushing airway obstructions deeper.

Universal Dual-Mask Sizing: Includes interchangeable adult and pediatric PVC face masks to serve diverse demographic markets.

Battery-Free Reliability: Purely mechanical suction design ensures zero maintenance and instant operation in high-stress rescue scenarios.

OEM/ODM Customization Ready: Designed for private labeling, custom color matching, and custom medical packaging for global B2B clients



Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

How to Use

I **Step 1: PLACE** – Position the mask securely over the victim's nose and mouth, pressing firmly on the mask edges to ensure a complete, airtight seal against the face and prevent air leakage.

Step 2: PRESS – Grip the PE handle and press downward to compress the airbag and exhaust internal air; the built-in one-way valve design ensures that downward pressure does not push foreign objects deeper into the airway.

Step 3: PULL – Maintain a firm, sealed position against the face without lifting the mask, and pull upward on the handle. The rapid rebound of the airbag generates powerful negative pressure to extract foreign objects.

Step 4: REPEAT – If the obstruction is not cleared on the first attempt, keep the mask sealed to the face and continuously repeat the press-and-pull cycle until the airway is cleared.



Materials & Technical Selection & Engineering

Face Mask (Medical-Grade PVC): Features a soft, skin-friendly texture that comfortably conforms to various adult and pediatric facial contours to guarantee an airtight seal. Its transparent properties allow rescuers to monitor the patient's lips and facial color in real time.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Airbag / Bellows (PE): Designed with superior elasticity and durability, resisting aging and wear through repeated compression cycles to reliably generate consistent suction power.

Colored End Component (ABS): Constructed from high-rigidity, impact-resistant ABS that resists shattering, features a striking color for rapid identification and handling, and ensures robust structural connections.

Packaging (PE Film): Individually sealed in moisture-proof, dust-proof PE film bags to maximize shelf life and hygiene for retail or bulk distribution.



Application Scenarios & Use Cases

Household Emergency Readiness: An essential safety device for family living rooms, kitchens, and bedrooms to protect household members of all ages from unexpected choking emergencies.

Automotive & Travel Safety: Compact and portable design ideal for storage in vehicle glove boxes, travel bags, or road trip emergency kits for on-the-go protection.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

School & Institutional Protection: A reliable, battery-free safety addition for classrooms, cafeterias, daycares, and campus health stations where quick, non-invasive airway clearance is critical.

Public Dining & Hospitality: Recommended for restaurants, hotel dining halls, and hospitality venues to support staff response readiness during dining incidents.

Professional & Commercial Distribution: Perfect for medical distributors, safety equipment wholesalers, and e-commerce brands seeking high-demand, OEM-ready household safety gear.

PROTECTING YOUR LIFE
— Essential for Every Scenario —

 Family

 School

 Travel

 Elderly Homes

Technical Specifications

Item	Specification
Set Specification	1 Set
Product Color	Green

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)

Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

Product Size (CM)	Main body: 189.89.8CM Adult mask: 12.510.07.5CM Child mask: 107.56CM
Package Size (CM)	25*25*10CM
Material	Handle: PE, Mask body: PVC, Air bag: PE, Color part: ABS
Net Weight	174g
Gross Weight (with package)	185.1g
Packaging Material	PE Bag
Carton Size (CM)	51*38*51
Carton Net Weight (KG)	9.45KG
Carton Gross Weight (KG)	10.5KG
Qty per Carton	50 Sets
MOQ for Custom Color	1500 Sets
MOQ for Custom Logo	1500 Sets
Invoice	Inform in advance: 1% ordinary invoice / 13% VAT invoice

Patented Choking Rescue Device Set – Dual-Tube Manual Suction System

Our production lines feature official Our products are backed by robust international patent protection and comprehensive global compliance certifications, including US design patents, EU appearance patents, and rigorous international standard testing.

Sourcing a patented dual-tube choking rescue device protects your supply chain from market infringement and ensures full legal compliance.

As an experienced OEM anti-choking device manufacturer, we combine independent patents with FDA compliance to deliver reliable emergency solutions.

Website: www.daqida.com **WhatsApp / WeChat:** [+86 187 5767 2042](https://wa.me/8618757672042)


Email: sean@daqida.com **Address:** Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province



DAQIDA | Yiwu Daqida Plastic Products Co., Ltd.

PRODUCT CERTIFICATIONS

Safety • Quality • Compliance

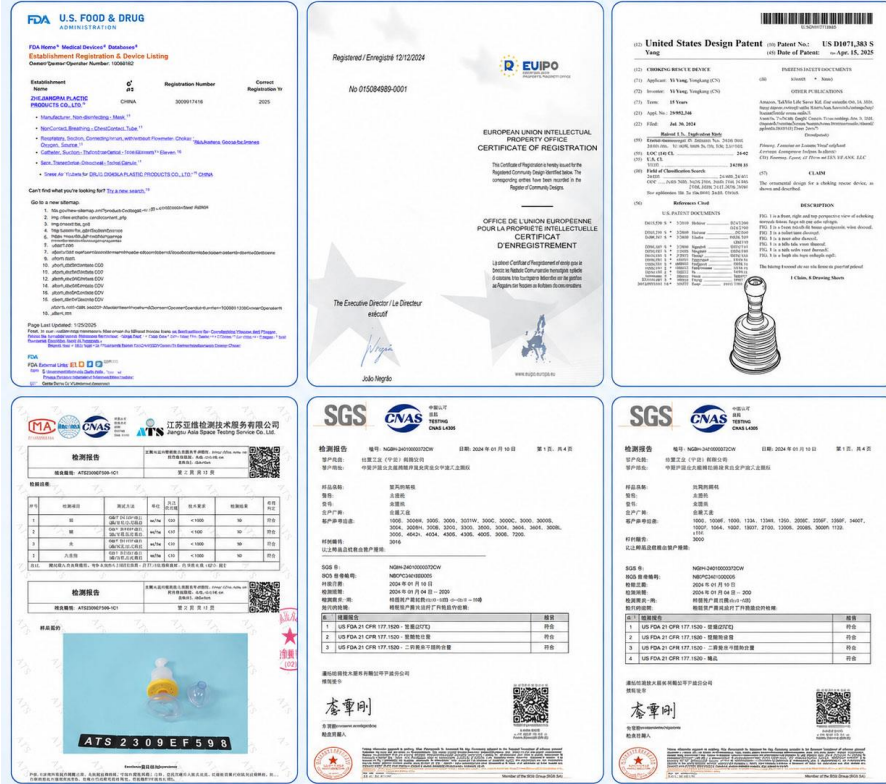

Website: www.daqida.com WhatsApp / WeChat: +86 187 5767 2042

Email: sean@daqida.com Address: Building 1, No. 47 Longxiang Road, Nansanli, Yiwu City, Zhejiang Province